




# WFCCD / WF4K

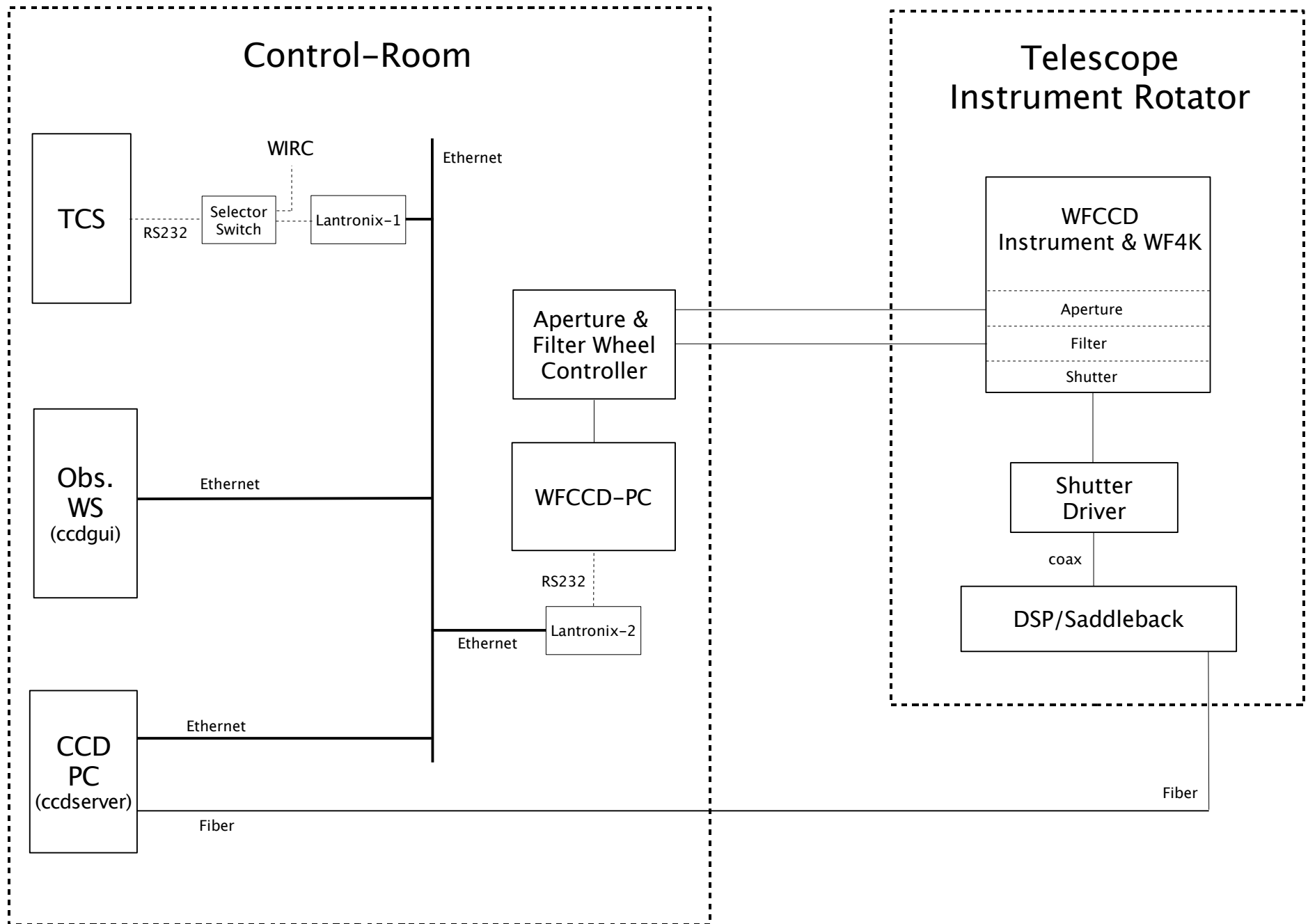


# Direct CCD / SItE2K

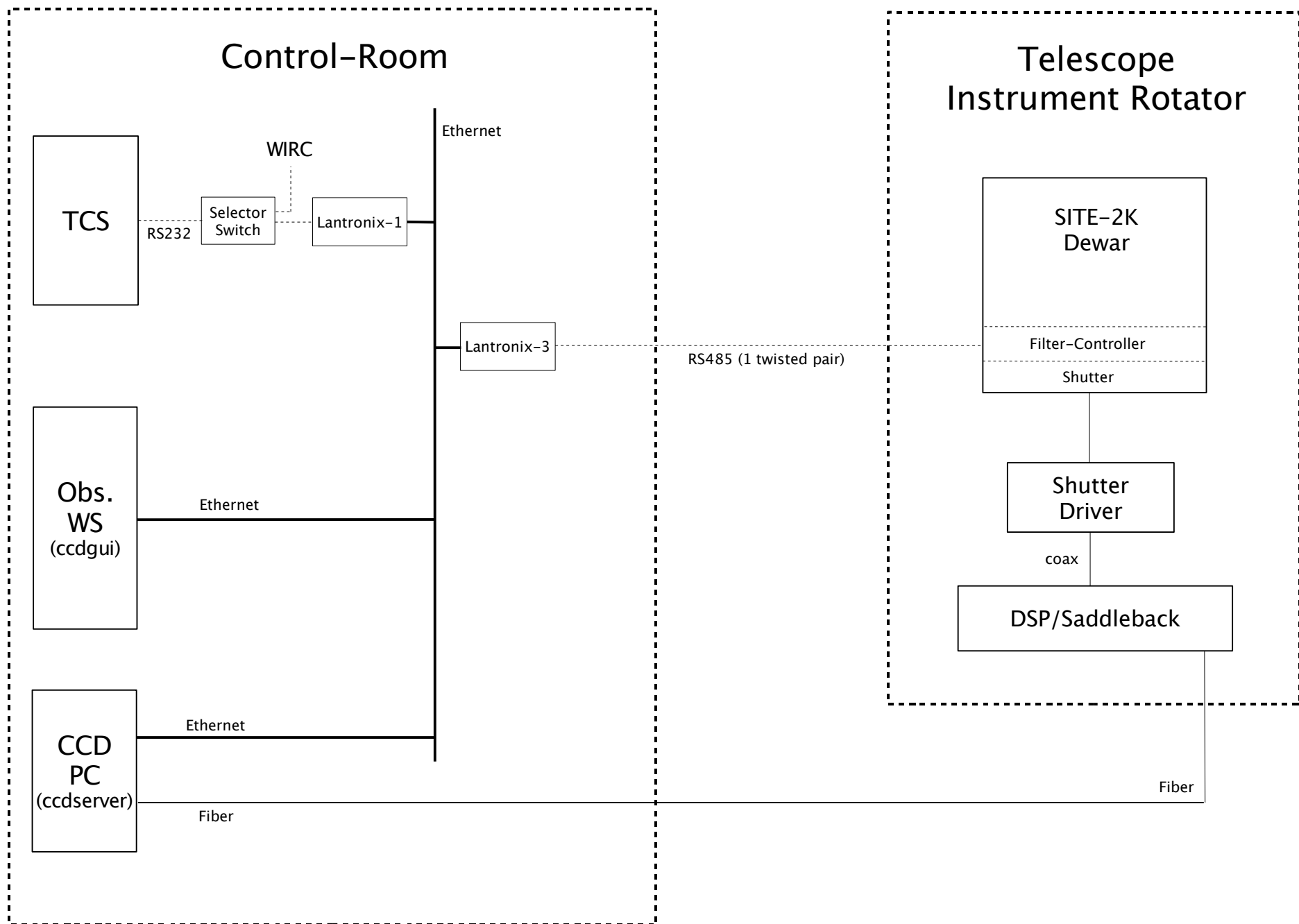


# Echelle / SITe2K

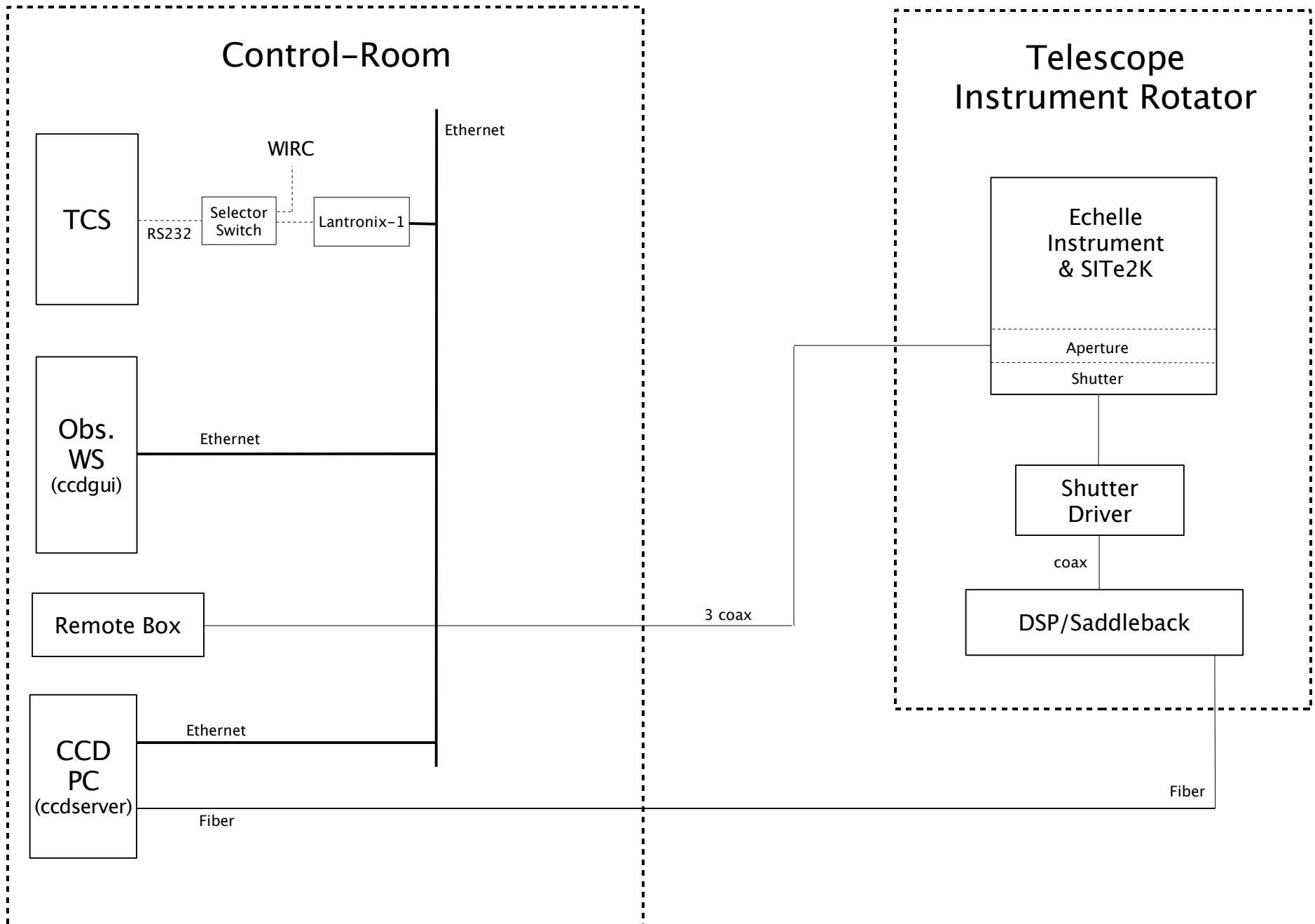

# B&C

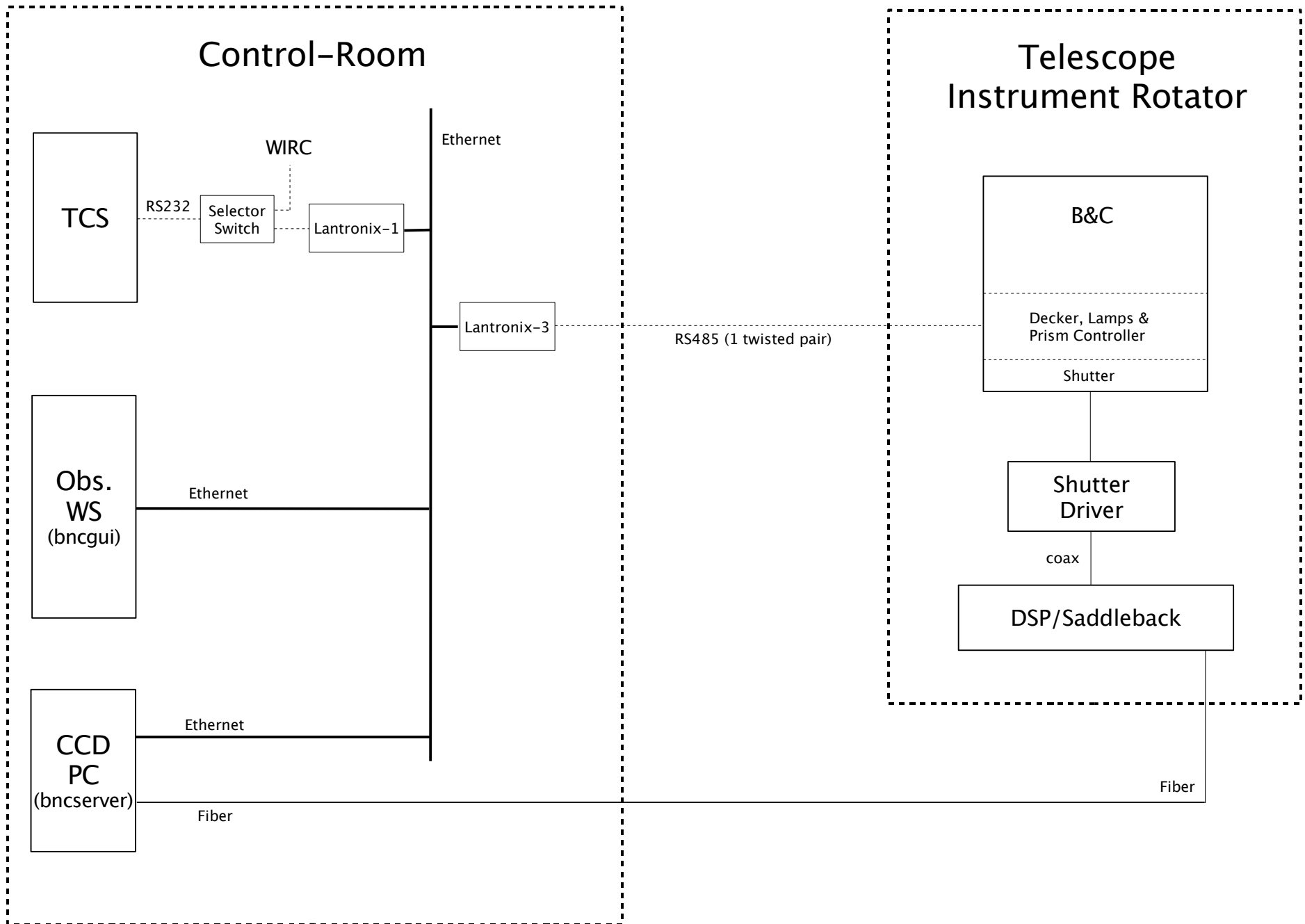

# CAPS-Cam

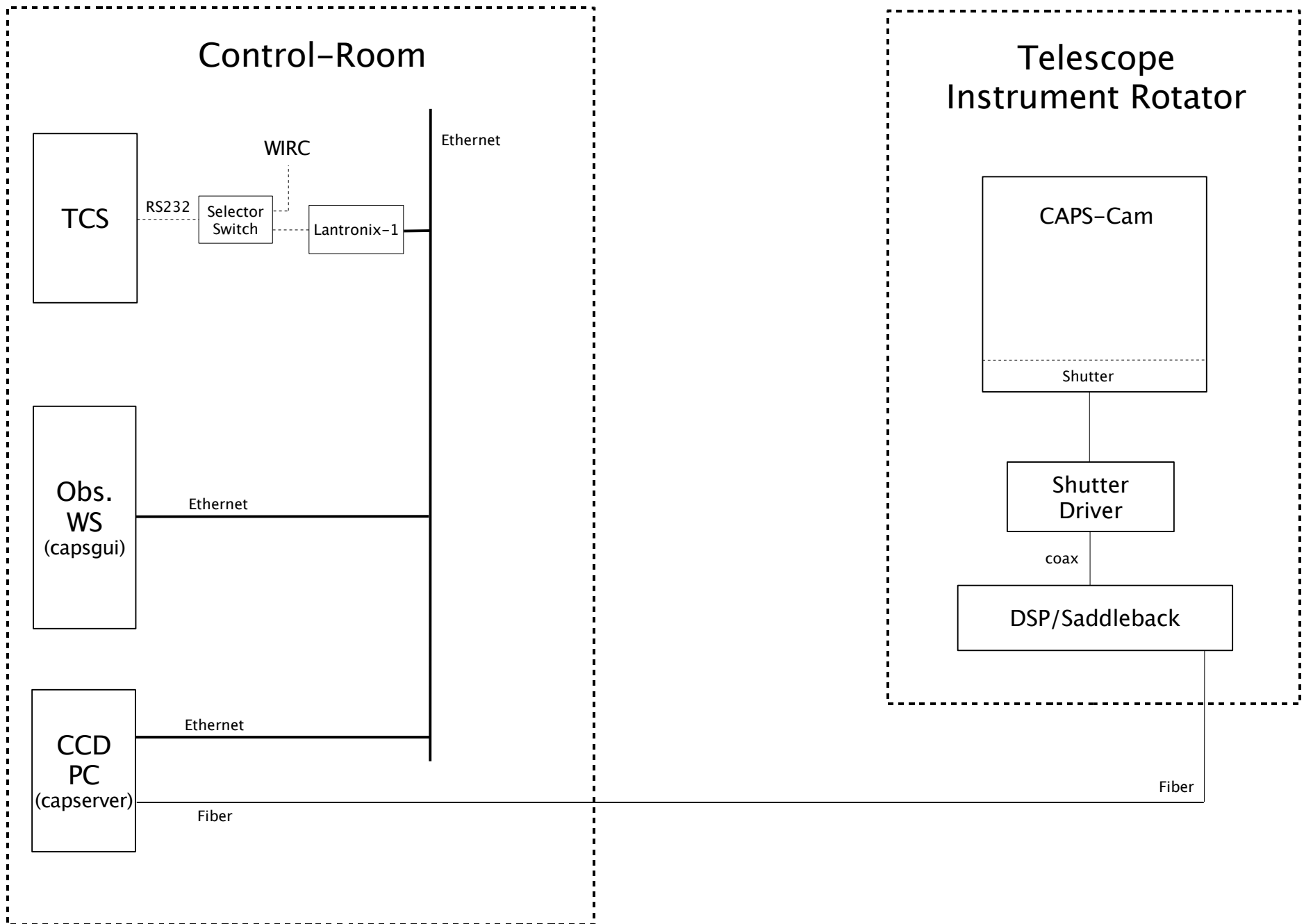

## Lantronix Setup

Lantronix	IP-Number	Parameters	Interface	Flow-Control	Connect	Disconnect	Flush
1 dupont-tcs	139.229.101.106	RS-232, 19200-8-1-N	0x4C	None	0xC0	0x00	0x00
2 dupont-wfccd	139.229.101.107	RS-232, 9600-8-1-N	0x4C	None	0xC0	0x00	0x00
3 dupont-imb	139.229.101.108	RS-485, 19200-8-1-N	0x4F	None	0xC0	0x00	0x00

## Software Changes

CCD: (WFCCD, Direct, Echelle):

- \* Magellan style 'ccdgui' for WFCCD/Direct/Echelle instruments running on the observer workstation
- \* Magellan style 'ccdserver' (SITE-2, WF4K chips) running on the CCD-PC
- \* TCS- communication via Lantronix (#1) TCP/IP-to-serial converter
- \* WFCCD-communication via Lantronix (#2) TCP/IP-to-serial converter
- \* Direct-Filter communication via Lantronix (#3) TCP/IP-to-serial converter

BNC:

- \* bncgui: TCS-communication via Lantronix (#1) TCP/IP-to-serial converter
- \* bncgui: to include the bnc.tcl script functions and use the Lantronix (#3) converter

CAPS:

- \* capsgui: TCS-communication via Lantronix (#1) TCP/IP-to-serial converter

## Hardware Requirements

- \* 3 Lantronix (UDS-10) TCP/IP-to-serial converters (already at SBS)
- \* Adapter cable to TCS: DB25<->DB9 (?)
- \* Adapter cable to WFCCD: DB25<->DB9 (?)
- \* Adapter cable to Direct-Filter/BNC: DB25<->DB9 (?)